



FEATHER SURGICAL BLADE

No.350 / No.391

No.350


The blade shape allows to perform a micro incision.

Width: 2.0mm

No.391

Allows to perform a safer incision and high operability.

Total length is 3mm longer and blade length is approximately a half of No.390.

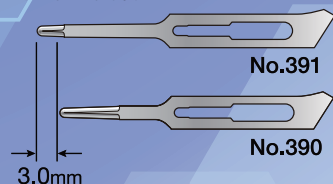

Feather surgical blades No.350/No.391 are the small sized surgical blades that can be handled in a narrow space of the oral cavity. These blades enable to perform the micro incision especially for the microsurgery. Please use No.350/No.391 as well as other Feather surgical blades No.370, No.390, and No.390C for the micro surgery according to the cases.

The blades can attach to Feather Surgical Handle.

Applicable handles

No.3



No.7




DF-160S-3




MADE IN JAPAN

Example use of Feather Surgical Blade No.350/No.391

No.350

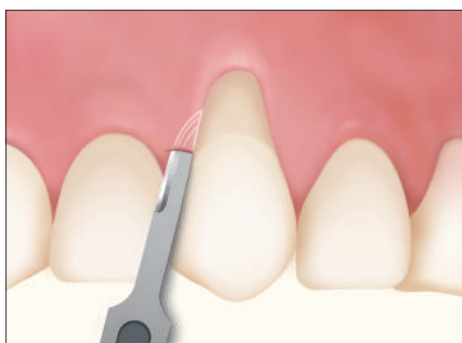


It allows the delicate incision in the thin gingiva and creates an accurate flap design.

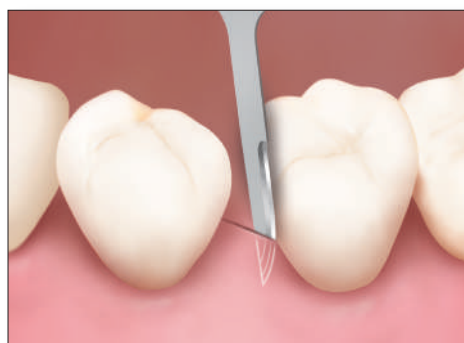


A small sized surgical blade is reachable between the teeth during the periodontal regeneration surgery; allows an oblique incision in the interdental space and creates the accurate flap.

No.391

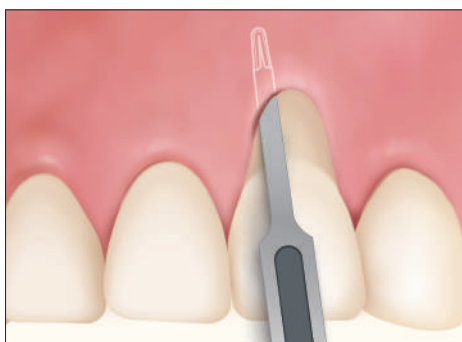

Total length is longer than No.390 and it allows a deeper envelope flap.

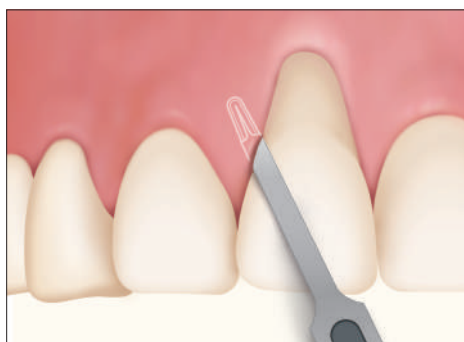

A shorter blade than No.390 allows a safe and accurate gingiva incision during the envelope flap surgery.

| Blade shape | Blade thickness | Cutting edge | Package |
|-------------|-----------------|---------------------|-------------------|
| No.350 | 0.38mm | Sharp on both sides | 10 blades per box |
| No.391 | | Sharp on one side | |
| No.370 | | Sharp on both sides | |
| No.390 | | Sharp on one side | |
| No.390C | | Sharp on both sides | |



FEATHER SAFETY RAZOR CO., LTD.
OVERSEAS TRADE DIVISION

3-70, OHYODO-MINAMI 3-CHOME, KITA-KU, OSAKA, 531-0075 JAPAN
PHONE: +81-6-6458-1638 FAX: +81-6-6458-1611
URL: www.feather.co.jp/en/ E-mail: overseas@feather.co.jp



ISO 13485 Certified