



SURGERY

- Surgical blade and handle
- Blade remover
- Blade pliers
- Disposable scalpel
- Safeshield scalpel
- Blood lancet
- Disposable dermatome
- Skin graft knife blade
- Biopsy punch
- Skin scraper
- Chokkin
- Disposable prep razor



MADE IN JAPAN

FEATHER SURGICAL BLADE

Feather surgical blades are available in either stainless steel or non-glare carbon steel that satisfies a wide range of demands of surgeons. The blades are individually packed and sterilized.

FEATHER is a leader in technological innovation and the creation of extremely sharp, durable blades.

	Stainless Steel	Carbon Steel
10		
11		
12		
12d		
14		
15		
15s		
15c		
20		
21		
22		
23		
24		
25		

Actual size

Packaging: Box of 100 pcs.

Stainless Steel		Actual size	
No.350		No.390	
No.370		No.390C	
		No.391	

Feather surgical blade No.350,370,390,390C and 391 blades are small shaped surgical blades, suitable for microsurgery (fine incision) with microscope. With appropriate shape and thickness, these blades are easy to control and to use for fine incision. Compatible with Feather Surgical Blade Handles (No.3, 7 and DF-160S-3).

Packaging: Box of 10 pcs.

FEATHER BLADE HANDLE FOR SURGICAL BLADE

Various types are available to meet diverse needs. Each is scientifically designed for smooth loading, handling ease and perfect balance. The handles are made of stainless steel (except DF-160S-3) to assure long life service and the blade integrates into the handle to form a single instrument. Handles No.3, 7 and DF-160S-3 accommodate blades No.10, 11, 12, 12d, 14, 15, 15c, 15s, 350, 370, 390, 390C and 391. Handle No.4 is designed for use with blades No.20, 21, 22, 23, 24 and 25.

No.3



No.7



DF-160S-3



Round-shaped handle made of titanium alloy. With smooth round shape, great weight balance and ideal length (160mm), this handle provide superior operability for incision in fine/deep area.

DLC coating

DF-160S-3 is coated with a diamond like carbon film (DLC) that offers superior performance and durability. The carbon coating prevents the handle from rusting and corroding and the handle design allows for smooth loading of blades.



No.4



FEATHER BLADE REMOVER

Feather Blade remover enables safe and easy removal of used disposable blades from any size or type of scalpel handles. It improves safety and hygiene by eliminating direct contact with blades and facilitates replacement.

Simple - - Three steps.

Safe & Hygienic - - No direct contact with blades for protection of staff against infection.

Disposable - - When this economical remover becomes full of used blades, discard it.



● Non-sterile

	Blade remover
Dimension(mm)	58(H) x 79.3(W) x 25(D)
Weight(g)	21.7
Blade insertion slot(mm)	17 x 5.2
Material	Polypropylene



- Sterile
- Individually packaged
- Just discard after use.

	Blade remover S
Dimension(mm)	19.3(H) x 79(W) x 27(D)
Weight(g)	11.8
Blade insertion slot(mm)	17 x 5.2
Material	Polypropylene

FEATHER BLADE PLIERS

Feather Blade Pliers No.3 enables users to remove a Surgical Blade from the Handle safely. The pliers can be cleaned and resterilized after use.



<CAUTIONS>

● Exclusively for Feather Surgical Blade Handle No.3, 7 and DF-160S-3. (Cannot be used for No.4.)

● Cannot be used for surgical blade handles of other manufacturers.

	Blade pliers
Weight(g)	40
Material	Stainless steel
Packaging	1 pc.

FEATHER DISPOSABLE SCALPEL

The Feather disposable scalpel is the last word in precision. Feather, utilizing feedback from surgeons, has developed a design with the sharpest cutting edge available today. Well-balanced construction, remarkable durability and in-the-package sterility further enhance total performance. Best of all, these throw-away scalpels are highly economical.

Economical & efficient

Fully automated Feather production and quality control systems are the reasons why these disposable scalpels are the first choice of surgeons everywhere. Feather disposable scalpels provide unsurpassed quality, uniformity and economy.

Indispensable in emergencies

Each scalpel is sterile and ready for use as soon as it is removed from the pull-apart aluminum wrapper. You can carry and use them.

Disposable scalpel



70% size

Disposable scalpel <Deep>

These scalpels feature a plastic handle about twice as long as the standard type for convenient use in difficult-to-reach areas. Available with blades No. 10, 11, 15, 20, 21, 22 and 23.



70% size

Disposable scalpel <Mini>


This type features a slender plastic handle and is highly suitable for use in surgery in limited areas. Available with blades No. 11, 14 and 15.



70% size


	Disposable scalpel	Disposable scalpel <Deep>	Disposable scalpel <Mini>
10			
11			
12			
12d			
14			
15			
15c			
20			
21			
22			
23			
24			
25			

How to open the package

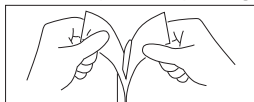

Simply peel the package off by pulling apart the two flaps, as indicated.




How to remove the snap-off cover

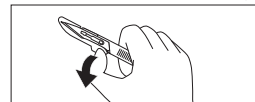



The blade cover is easily removed by applying firm pressure to the snap-tab, as indicated.

Safeshield Scalpel

Clear protective shield can be moved forward (to cover the blade) and retracted multiple times. 11 sizes are available and color-coded button on the protective shield assists blade identification.

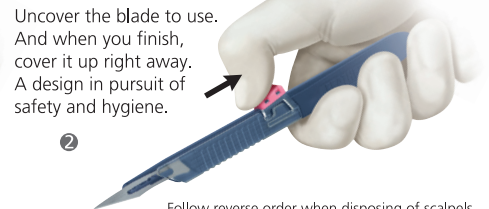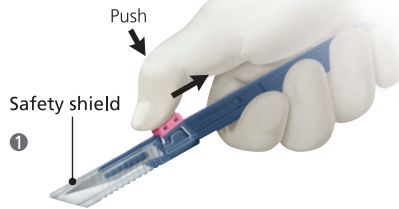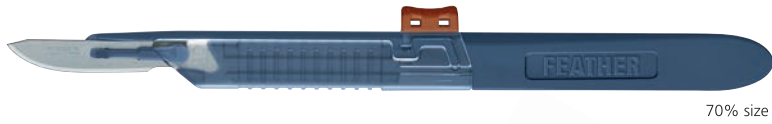

Follow reverse order when disposing of scalpels.

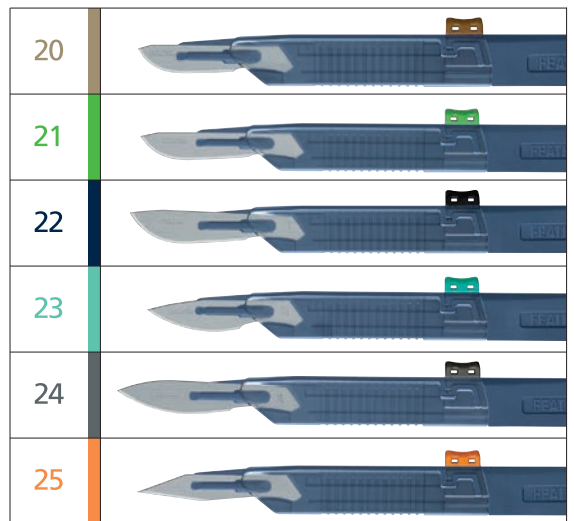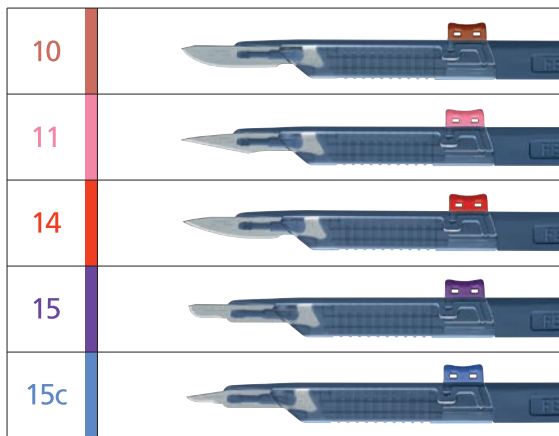



How to open the package



Simply peel the package off by pulling apart the two flaps, as indicated.


Safeshield Scalpel Vessel

Feather "Safeshield Scalpel Vessel" is designed for custom procedure tray / kit. The safety cover protects the blade in procedure and handling. Supplied in plastic box of 100 pcs, non-sterile.

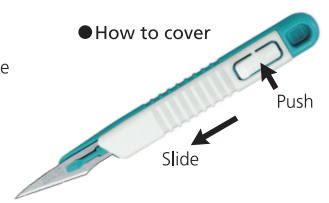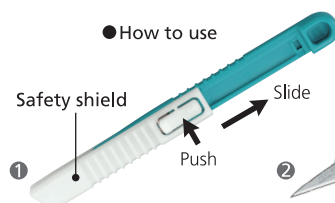

FEATHER BLOOD LANCET


Feather Blood Lancet assures superior hardness and sharpness. Each lancet is individually packed and sterilized.

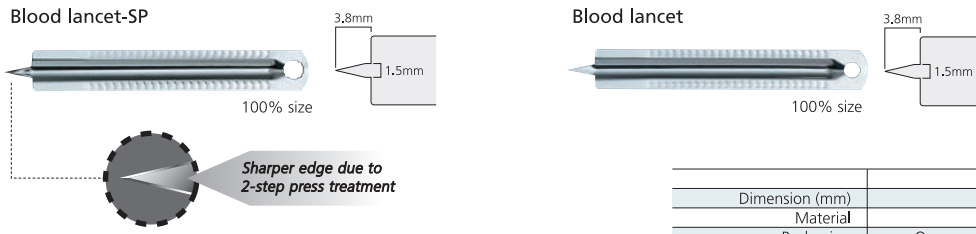

	Blood lancet-SP/lancet
Dimension (mm)	40(L) x 5.4(W)
Material	Stainless steel
Packaging	One small box contains 200 pcs.

DERMATOLOGY

FEATHER DISPOSABLE DERMATOME

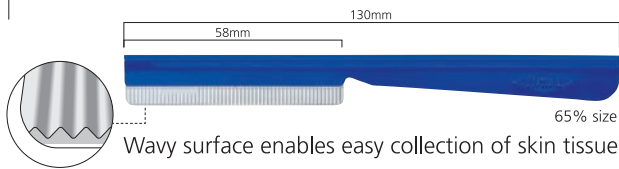

		Disposable dermatome
Dimension	Blade (mm)	58(L) x 16(W)
	Handle (mm)	130(L) x 12(W)
Material	Blade : Stainless steel (specially processed blade)	
Packaging	Individually packed and sterilized. Box of 20 pcs.	

FEATHER SKIN GRAFT KNIFE BLADE


Feather skin graft knife blade is a disposable blade made of stainless steel. Its sharp edge offers a smooth grafting for the burnt skin. Skin graft knife blade (SGB-157) is compatible with the following knives: ROSENBERG, BRAITHWAITE, COBBETT, and WATSON.



		Skin graft knife blade(SGB-157)
Dimension (mm)	157(L) x 18(W) x 0.229(T)	
Material	Stainless steel	
Packaging	Individually packed and sterilized. Box of 20pcs.	

FEATHER BIOPSY PUNCH

Feather Biopsy Punch offers excellent sharpness by ultra-precision sharpening technology which originates in surgical blade manufacturing.

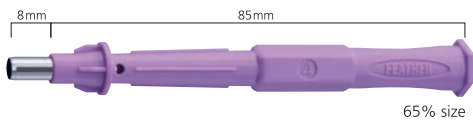

		Biopsy punch
Inner Diameter (mm) / Product No.	1.5 / FBP-15	3.0 / FBP-30
	4.0 / FBP-40	5.0 / FBP-50
Material	Blade	
	Stainless steel	
Packaging	Handle	
	Plastic	
	Individually packed and sterilized. Box of 20 pcs.	

FEATHER SKIN SCRAPER

Feather Skin Scraper is a spatula suitable for tissue collection for diagnosis of mycosis. Edge is not sharpened and it can be safely used without cutting skin. 1mm thick spatula allows firm and easy tissue collection. Autoclavable.



	Skin Scraper
Length(mm)	140
Material	Stainless steel
Packaging	1pc.

Chokkin (with strap)

Pocket-sized scissors designed for exceptional handling ease, cutting ability and safety. "Chokkin" scissors are small enough to fit in a nurse's pocket. And come with a holder and clip for easy access without damage to the pocket. The perfect tool for the day-to-day hospital routine of nurses. Blade edges are specially processed with fluoroplastic. Residual glue from adhesive tape will not cling to blade after cutting, thus assuring smooth cutting reliability of the scissors.

Improved Safety


The blade tip is rounded and protected with a plastic guard. When cutting bandages, only the plastic guard touches the skin. Thus, adding relief and reassurance to the patient during cutting procedures.

Friendly Design

"Chokkin" are designed with gentle, curved lines.



Elastic strap provided

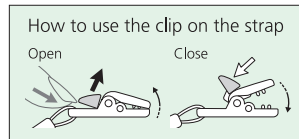

*Stainless steel blade with plastic handles (ABS resin, heatproof temperature : 70°C).

Disposable prep razor 121S-B

*Quick and easy to use. *Individually blister packed. *Superbly sharp double-edge blade.



		Disposable prep razor 121S-B
Handle length(mm)		65
Material	Blade	Stainless steel
	Handle	Plastic
Packaging		Supplied in box of 50 pcs.



FEATHER SAFETY RAZOR CO., LTD.
OVERSEAS TRADE DIVISION

3-70, OHYODO-MINAMI 3-CHOME, KITA-KU, OSAKA, 531-0075 JAPAN
PHONE: +81-6-6458-1638 FAX: +81-6-6458-1611
URL: www.feather.co.jp/en/ E-mail: overseas@feather.co.jp



ISO 13485 Certified

